



Technical Data Sheet

Diagenode sa

CHU, Tour GIGA B34, 3^e étage
Avenue de l'Hôpital, 1
4000 Liège - Belgium

Product name:
antibody directed against hLeo1

(Component of the human PAF complex)

Catalog #: pAb-017-050	Type: Polyclonal	Size: 50 µg/ 25 µl
Lot #: 001	Source: Rabbit	Concentration: 2.0 µg/µl

Description: This antibody has been raised against the full length recombinant human Leo1 protein.
PAF is a five-subunit protein complex composed of Paf1, Cdc73, **Leo1**, Rtf1 and Ctr9, which was first purified from yeast in association with RNA polymerase II and which is believed to function in transcription elongation [1]. See overview below.

Specificity: Human: positive
Other species: not tested

Applications	Suggested dilution	References
ELISA	Not tested	
Dot blotting	Not tested	
Western blotting	1:2000	Fig 1
Gel Supershift	Not tested	
Immunochemistry	Not tested	
Flow cytometry	Not tested	
Immunoprecipitation	Not tested	
ChIP	1:100	[2]

Format: In solution in PBS including 0.05% azide and 0.05% ProClin 300. The polyclonal antibody has been protein G purified.

Storage: For long storage, store at -20°C. Do not freeze-thaw.

Precautions: This product is for research use only. Not for use in diagnostic or therapeutic procedures.

References:


[1]. Rondon A.G., Gallardo M., Garcia-Rubio M. and Aguilera A. 2004 *EMBO Rep.* 5(1):47-53.

[2]. Zhu B., Mandal S.S., Pham A.D., Zheng Y., Erdjument-Bromage H., Batra S.K., Tempst P. and Reinberg D. 2005 *Genes Dev.* 19:1668-73.

Availability date: March 12, 2007. Last data sheet update: March 19, 2007

Lot #: 001: rabbit#1/ bleed day: final bleed/ purification day: February 28, 2007

Figure 1:

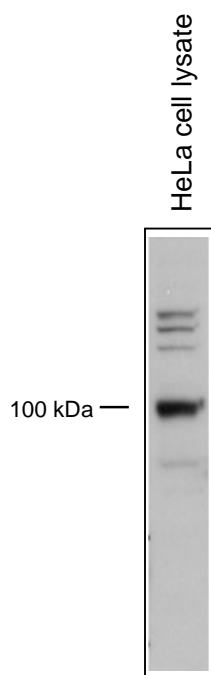

Western blot analysis using the antibody from Diagenode anti-hLeo1

HeLa cell lysate was analysed by Western blot using the antibody directed against hLeo1 (Diagenode, crude serum was used for this WB shown here above), used at dilution 1:2000. Size and location of the protein are indicated.

Overview

The yeast PAF (yPAF) complex interacts with RNA polymerase II and coordinates the setting of histone marks associated with active transcription. The isolation and functional characterization of the human PAF (hPAF) complex was reported. hPAF shares four subunits with yPAF (hCtr9, hPaf1, hLeo1, and hCdc73), but contains a novel higher eukaryotic-specific subunit, hSki8 [2]. In addition to coordinating events during transcription (initiation, promoter clearance, and elongation), hPAF also coordinates events in RNA quality control [2].